

Safe2 Rotator™ Luer Connector


The connection you can rely on

Easy attach – Firm lock – 360° Rotational flexibility - Easy detach

Good and safe connection of medical tubing is critical to safe and effective therapy. Tube leakages and unstable connections can result in serious and irreversible patient injury or even death.

When it comes to stopcocks, apart from being securely connected, they must be convenient to use and perform their primary functions – direct fluids and provide access to the administration set.

With this in mind, Elcam Medical has designed the new Safe2 Rotator™.

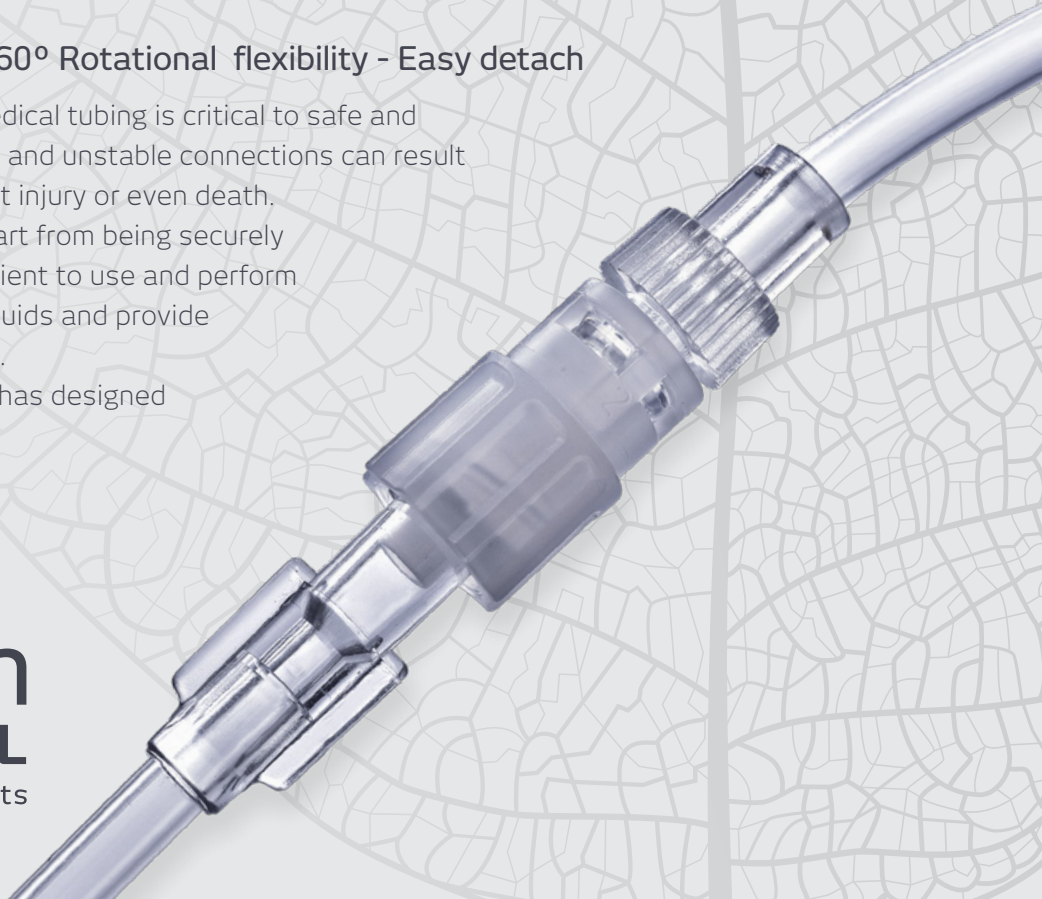

Secured Connections



 Elcam Devices


 Elcam Stopcocks

 Elcam Accessories

Safe2 Rotator™ Luer Connector

The connection you can rely on

The Safe2 Rotator™ flexible material (polyethylene) and spinning lock feature provide ease of use and an extra layer of passive protection to your most critical fluid managing applications, without altering the traditional usability and common practice:

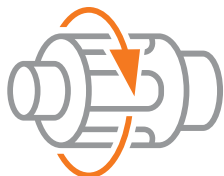


Mechanical advantages

- Reduces mechanical stresses induced during prolonged connection and minimizes crack formation and consequent leakages
- Prevents leakages due to over-tightening
- Prevents "stuck" luer locks due to adhesion and facilitates easy disconnection upon end of therapy
- Prevents accidental disconnections



Safety advantages

- Minimizes drug and blood leakages and related adverse events:
 - Blood borne pathogens exposure
 - Hazardous drug exposure
 - Blood lose and dehydration
 - Air embolism
 - Microbial ingress
 - CRBSI
- Maximizes therapy completeness

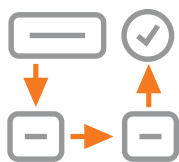

Workflow advantages

- Reduces tube kinks and twists in the tubing set
- Facilitates convenient approach to injection and sampling ports during fluid administration
- Enables rotation of coupled luer lock connections and stopcocks

Available configurations of the connector (stand-alone version):

Tube OD (mm)	Body Material	NUT Material	NUT color
6.8	PC	PE	○ ● ● ● ●
4.1	PC	PE	○ ● ● ● ● ●



Especially convenient in use with needleless connectors providing rotational flexibility of the distal part of the set.

Elcam Medical's Safe2 Rotator™ – secured connection along the entire treatment

The Safe2 Rotator™ Luer Connector is available in both stand-alone and stopcock configurations. Stopcock configuration is CE approved for non-sterile products and 510(k) cleared

[Click here to contact us >>](#)



www.elcam-medical.com

Head Office/Manufacturer
Elcam Medical

BarAm 1386000,
Israel
Tel: 972 4 698 8120
Fax: 972 4 698 0777
sales@elcam.co.il

USA
Elcam Medical Inc.

1 University Plaza, Ste. 511
Hackensack, NJ 07601
Tel: +1 201 4571120
Fax: +1 201 4571125
info@elcam-medical.com

Europe
Elcam Medical Europe

Via Emilia Romagna 15,
41012 Carpi (Modena) Italy
Tel: +39 059 697710
Fax: +39 059 697705
sales.eu@elcam-medical.com

China
Elcam China Office

Room 3309, Jinhai Mansion
No. 89 of Anshun Road,
Shanghai, PRC. Zip: 200051
Tel: +86 21 62829967
Fax: +86 21 62829967
johnliu02@126.com